

MAA573Hu22

Monoclonal Antibody to Procollagen III N-Terminal Propeptide (PIIINP)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)



## [ PROPERTIES ]

**Source:** Monoclonal antibody preparation

Host: Mouse

Antibody isotype: IgG2b Kappa

**Purification:** Protein A + Protein G affinity chromatography

Clone number: E10

Traits: Liquid

Concentration: 1mg/mL

**UOM**: 200µg

Applications: WB; IHC; ICC; IP.

#### [ IMMUNOGEN ]

Immunogen:

**Accession No.:** 

# [ APPLICATIONS ]

Western blotting: 0.5-5µg/mL

Immunohistochemistry: 5-30µg/mL

Immunocytochemistry: 5-30µg/mL

Optimal working dilutions must be determined by end user.

#### [FORMULATION]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN3, 50% glycerol.

## [QUALITY CONTROL]

**Content:** The quality control contains recombinant PIIINP disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL).

**Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.



**Loading Buffer:** 100mM Tris (pH6.8), 1% SDS, 150mM NaCl, 50% glycerol, 0.02% BPB, 50mM DTT and 0.02% NaN $_3$ .

#### [STORAGE AND STABILITY]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

### [ IDENTIFICATION ]

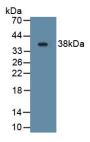


Figure. Western Blot; Sample:

Recombinant PIIINP, Human.

### [ IMPORTANT NOTE ]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.