

PAA573Hu01

Polyclonal Antibody to Procollagen III N-Terminal Propeptide (PIIINP)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)



# [PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

**Purification:** Antigen-specific affinity chromatography followed by Protein A affinity

chromatography

Traits: Liquid

Concentration: 500µg/mL

**UOM**: 200µg

Applications: WB; ICC; IHC-P; IHC-F; ELISA.

[ IMMUNOGEN ]

Immunogen:

**Accession No.:** 

### [APPLICATIONS]

Western blotting: 0.5-2µg/mL;1:500-2000

Immunohistochemistry: 5-20µg/mL;1:50-200

Immunocytochemistry: 5-20µg/mL;1:50-200

Optimal working dilutions must be determined by end user.

#### [FORMULATION]

**Form & Buffer:** Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

## [QUALITY CONTROL]

**Content:** The quality control contains recombinant PIIINP disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL).

**Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris (pH6.8), 1% SDS, 150mM NaCl, 50% glycerol, 0.02% BPB,

50mM DTT and 0.02% NaN<sub>3</sub>.

#### [STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]



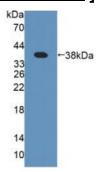


Figure. Western Blot; Sample: Recombinant PIIINP, Human.

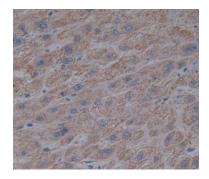

Used in DAB staining on fromalin fixed paraffin- embedded liver tissue

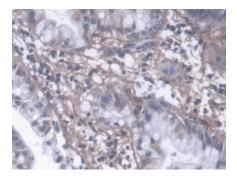

DAB staining on IHC-P;
Samples: Human Liver cancer Tissue

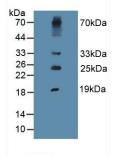


Figure. Western Blot; Sample: Porcine Liver Tissue.

# [ IMPORTANT NOTE ]

The kit is designed for research use only, we will not be responsible for any issue if the kit was



used in clinical diagnostic or any other procedures.