

RPB316Hu02 100µg

Recombinant Procollagen C Proteinase Enhancer 1 (PCPE1)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)



### [PROPERTIES]

**Source:** Prokaryotic expression

Host: E.coli

Residues: Gly38~Asp449

Tags: N-terminal His Tag

Subcellular Location: Secreted

**Purity:** > 95%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% SKL, 5% Trehalose

and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.5

Predicted Molecular Mass: 47.8kDa

Accurate Molecular Mass: 48kDa as determined by SDS-PAGE reducing conditions.

#### [USAGE]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

#### [STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.


**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.



## [SEQUENCE]

			GGD	VKGESGYVAS
EGFPNLYPPN	KECIWTITVP	EGQTVSLSFR	VFDLELHPAC	RYDALEVFAG
SGTSGQRLGR	FCGTFRPAPL	VAPGNQVTLR	MTTDEGTGGR	GFLLWYSGRA
TSGTEHQFCG	GRLEKAQGTL	TTPNWPESDY	PPGISCSWHI	IAPPDQVIAL
TFEKFDLEPD	TYCRYDSVSV	FNGAVSDDSR	RLGKFCGDAV	PGSISSEGNE
LLVQFVSDLS	VTADGFSASY	KTLPRGTAKE	GQGPGPKRGT	EPKVKLPPKS
<b>QPPEKTEESP</b>	SAPDAPTCPK	QCRRTGTLQS	NFCASSLVVT	ATVKSMVREP
GEGLAVTVSL	IGAYKTGGLD	LPSPPTGASL	KFYVPCKQCP	PMKKGVSYLL
MGQVEENRGP	VLPPESFVVL	HRPNQDQILT	NLSKRKCPSQ	PVRAAASQD

## [ IDENTIFICATION ]

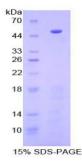

# [ IMPORTANT NOTE ]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.