

1kb DNA Ladder

Product Name	Product Code	Kit Packing**
1kb DNA Ladder	MBT051-50LN	50 Lanes (250 µl)
1kb DNA Ladder	MBT051-200LN	200 Lanes (4 x 250 µl)

**** The product is supplied with a vial of 6X Gel loading Buffer.**


Description:

The 1kb DNA ladder contains ten double-stranded DNA fragments ranging from 200bp to 10kb. The amounts of 500, 2000 and 5000 bp fragments are adjusted so that their increased intensities upon staining serve as easy reference after gel electrophoresis. All fragments are precisely quantified and mixed during the production. The ladder is pre-mixed with loading dye and is ready to use.

Recommended load: 5 µl / lane

Storage conditions:

The ladder should be stored at -20°C. When stored under the recommended conditions, the product is stable for 2 years.

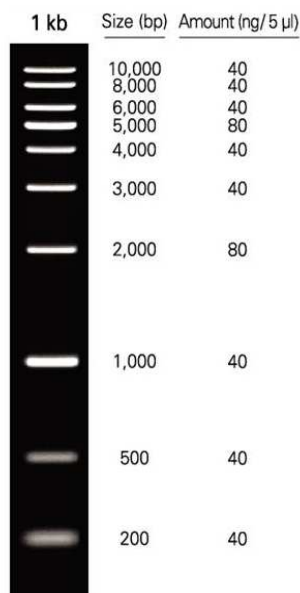

Quality Control:

Each lot of HiMedia's 1kb DNA Ladder is assayed for contaminating endonuclease, exonuclease, and non-specific DNase activities.

Technical Assistance

At HiMedia, we pride ourselves on the quality and availability of our technical support. For any kind of technical assistance, mail to mb@himedialabs.com

PIMBT051_0/0514

MBT051-04

Disclaimer :

User must ensure of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ Publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. Himedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal diagnostic or therapeutic use but for laboratory, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

CE A-516, Swastik Disha Business Park, Via Vadhani Indl. Est., LBS Marg, Mumbai - 400 086, India Customer Care No.: 022-6147 1919 Email: techhelp@himedialabs.com



Registered Office :

23, Vadhani Industrial Estate, LBS Marg,
Mumbai - 400 086, India.
Tel. : (022) 4017 9797 / 2500 1607
Fax : (022) 2500 2286

Commercial Office

A-516, Swastik Disha Business Park,
Via Vadhani Indl. Est., LBS Marg,
Mumbai - 400 086, India
Tel: 00-91-22-6147 1919
Fax: 6147 1920, 2500 5764
Email : info@himedialabs.com
Web : www.himedialabs.com

The information contained herein is believed to be accurate and complete. However no warranty or guarantee whatsoever is made or is to be implied with respect to such information or with respect to any product, method or apparatus referred to herein